

The *Clearview Pregnancy Test* has been chosen by the Sexual Health Resource Centre because we have found it is the most reliable, simplest to use, and the results are the easiest to interpret.

Please call us with any questions you might have in using the test.

Please remember that Queen's University Student Health Services (Queen's students only) or your family physician will administer a pregnancy test free of charge if you have OHIP, UHIP, or a similar health plan.

Please read through this pamphlet carefully before use.


Before Using the Test

- Do **NOT** use this test if you are taking the following drugs, as they will affect the accuracy of the test:
 - Profasi, APL, Follutin, PregnylNo other medications will interfere with this test.
- You must wait **14 DAYS AFTER** the last time you had unprotected sex or the date you suspect becoming pregnant in order to ensure accurate results.
- You must **USE YOUR FIRST MORNING URINE**.
- Perform this test **WITHIN 4 HOURS** of urine collection. If this is not possible, refrigerate the urine sample and then warm it to room temperature for 30 minutes before performing the test.

Instructions

1. Collect your **FIRST MORNING URINE** in the cup provided.
2. Take the pregnancy test out of its packaging.

3. Use the dropper provided to collect a small amount of urine from the cup. The dropper will fill with enough urine (at least 5 drops) for the test.
4. Make sure the pregnancy test is lying on a flat surface.
5. Carefully place at least 5 drops or more of urine into the sample well (see diagram). It is very important that you all the test to **LIE ON A FLAT SURFACE** while waiting to read the results. Handling the test could result in an inaccurate reading.
6. Wait 5 minutes, or until the control region shows a pink line. If the control region does not show a pink line, add another drop or two of urine to the sample well and wait another 5 minutes. If there is still no pink line in the control region, then there is a problem with this particular test. If this occurs, please contact us or come into the office and we will replace the test free of charge.

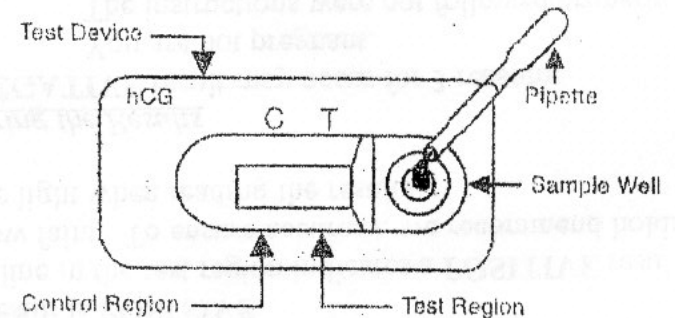

Reading Results

- If a pink line appears in the control region **ONLY**, the result is **NEGATIVE**.
- If pink lines appear in **BOTH** the control region and the test region, the result is **POSITIVE**.

- Any line in the test region indicates a **POSITIVE** result, regardless of how faint. To ensure accuracy, we recommend holding the test up to the light when reading the results.

Considering the Results

- A **NEGATIVE** result may occur for 3 reasons:
 1. You are not pregnant.
 2. The instructions were not followed properly. Repeat the test if you are concerned the test was performed incorrectly.
 3. You have not waited 14 days since you suspected becoming pregnant before performing the test. It is only after this period of time that the test can yield accurate results from your urine sample.
- A **POSITIVE** result means that you are pregnant. When the result is positive, the accuracy of the test is 100%. There are no false positives.
- If you test positive and then have your period, the test was not inaccurate; you most likely experienced a natural termination of your pregnancy. Up to one-third of all pregnancies end in spontaneous abortion and the bleeding caused by this natural termination appears very similar to regular menstrual flow. However, a small number of women still experience bleeding while pregnant. Therefore, if you experience bleeding or have a very weak positive result, we recommend that you confirm your results with another test after a few days or see a physician.

For more information about pregnancy options, support services available and birth control call the Sexual Health Resource Center at 533-2959.

**Sexual Health Resource Centre
533-2959**



**Pregnancy Test Information &
Instructions**